


## Packing Capillary Columns and Pre-columns (traps)

### PPE (personal protective Equipment)

Do not use any of the equipment without the appropriate training. Always wear appropriate PPE when working in the lab, including goggles, lab coats and gloves. PPE is provided, ask if you cannot find the appropriate PPE or if PPE is missing. Always wear goggles when working with fused silica. Review MSDS and SOP's before working with chemicals. Before working with the gas cylinders make sure you complete the EH&S "Compressed Gas Safety" online training (<https://www.ehs.washington.edu/training/compressed-gas-safety-online>).

### Packing Station Overview

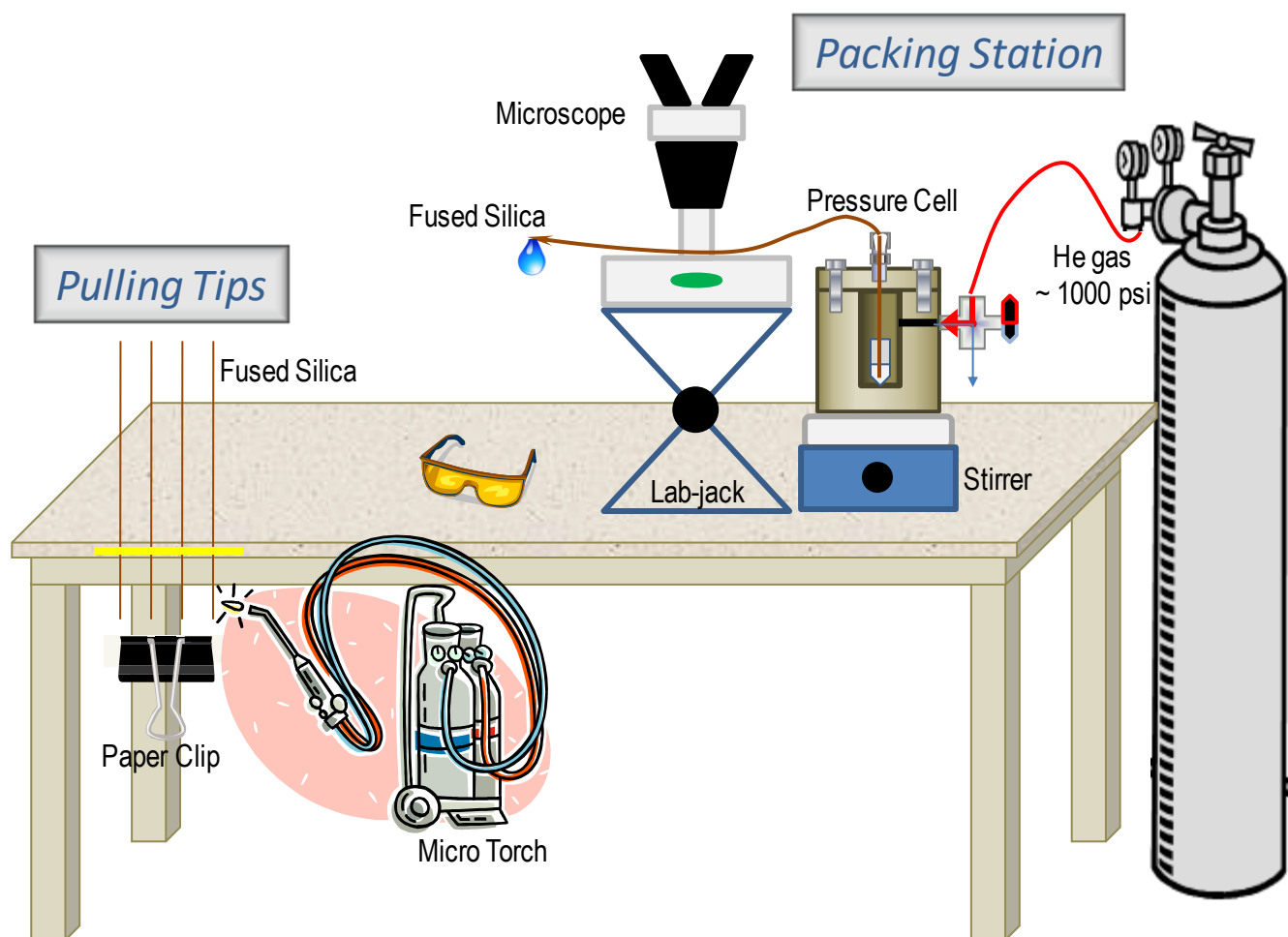

### Typical setup of a column packing station

The magnetic stirrer and the microscope are optional.

The table at the end of this document list parts needed to setup a packing station as well as some of the consumables. Note most parts can be obtained from multiple vendors, the table only lists examples.

The following sections describe how to pull tapered tips, make a frit and pack columns and traps.

For typical LC-MS setup see the LC-plumbing document ([https://proteomicsresource.washington.edu/docs/protocols05/LC\\_plumbing.pdf](https://proteomicsresource.washington.edu/docs/protocols05/LC_plumbing.pdf))

## Column/trap End Style

### Pulled tip (manual)

Depending on the resources available to you, you can use a laser tip puller or a micro-torch to pull a thin tip.


#### Reagents and Materials (see Table 1 below)

- Fused silica tubing 360  $\mu\text{m}$  OD by 50, 75, 100, 150, 200  $\mu\text{m}$  ID (Polymicro Technologies, see Table 1 below)
- Fused silica cutter
- Micro-torch (e.g. small torch for oxygen/propane (Kingsley North Inc.) or Microflame torch # 14606 (Alltech))
- OR Laser Fiber Puller (Sutter Instruments)

#### Procedure

Cut desired length of fused silica (5-10 cm longer than final column length) and use tape to attach the fused silica to the side of a table. Attach a large binder clip at the lower end of the fused silica. Use the micro torch to heat fused silica until it melts. The binder clip will act as a weight and pull the tip. Use a fused silica cutter to cut the tip to desired length and desired tip ID.

| SilicaTip™ Size vs. Flow Rate<br>(Recommendation from the NewObjective website) |                     |
|---|---------------------|
| Tip Size, ID ( $\mu\text{m}$ )  | Flow Rate* (nL/min) |
| 5   | 20-100              |
| 8   | 50-300              |
| 10  | 100-400             |
| 15  | 150-400             |
| 30  | 300-1000            |

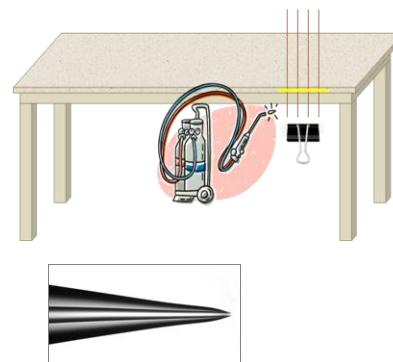

(Image from New Objective's Website)

### Pulled tip (Laser puller)

Alternatively, you can use a laser puller to pull your tip.

#### Reagents and Materials (see Table 1 below)

- Fused silica tubing 360  $\mu\text{m}$  OD by 50, 75, 100, 150, 200  $\mu\text{m}$  ID (Polymicro Technologies, see Table 1 below)
- Fused silica cutter
- Lighter, preferable a wind resistant lighter
- Laser Fiber Puller (Sutter Instruments P-2000)



Sutter instruments P-2000

**For First time users: Please ask for help when using the laser puller for the first time!**

#### Laser Puller Instructions for pulling $\leq 375 \mu\text{m}$ OD fused silica tips

1. Turn Power switch on left side of cabinet ON.  
To assure the most reproducible results you should allow the unit to warm u with the lid down for **15 min** before pulling tips
2. Type Program # and hit <Enter> key. On the UWPR instrument, my favorite is #1.

3. **Prepare fused silica:** the plastic sheath **MUST** be stripped to expose the quartz glass before loading into the puller. Use the flame of a lighter to burn off about 2 cm of the coating. Wipe CLEAN with a methanol soaked kimwipe.
4. Load the fused silica into the puller:
5. Close lid, be careful not to bread the fused silica
6. Press <Pull>. The glass should separate in a few seconds.
7. Loosen the clamping knobs to remove the fused silica

NOTE: if the fused silica is not aligned properly it will not separate and the red light will stay on. Press <Stop>, open the lid and re-align the fused silica. Close the lid and try again.

#### To select a different program

1. Hit the <Reset> button, wait for welcome screen to appear
2. Type program # and hit <Enter>

#### To create a new program

1. **Be courteous to others and do NOT overwrite someone else's programs!**
2. Select an unused program (see spreadsheet)
3. enter the values as desired.
4. hit <Pull>
5. If you wish to keep your program, indicate so on the spreadsheet by writing down the values of your program

#### HEAT (laser power output):

Range ~200-350 (if your program and glass require HEAT greater than 350 to get separation, then there is a problem. Ask for Help!  
Generally, higher HEAT tends to give longer finer tips

#### FIL (Filament, scan length):

Range 0-5; This is the scan length of the laser. The manufacturer recommends to use 0 for fused silica. However heating up a longer stretch of fused silica helps generating longer tips. Note extending the scan length may require higher HEAT to get the glass to separate.

| Filament # | Scan Length |
|------------|-------------|
| 0          | 1 mm        |
| 1          | 1.5 mm      |
| 2          | 1.9 mm      |
| 3          | 4.5 mm      |
| 4          | 6.5 mm      |
| 5          | 8 mm        |

#### VEL (Velocity, trip point):

Range 15-35; Changing the Velocity will affect the thickness of the glass. Generally if the glass is thinned too much, reduce heat (increments of 10) and increase the velocity (increments of 2), and vice versa.

#### DEL (Delay)

Range 0-255 (>128 recommended); Controls the timing of the start of the hard PULL relative to the h deactivation of the laser. If the delay is set to 128 the hard pull is initiated at the same time as the deactivation of the laser. If PULL set to 0, the delay allows the glass to cool before the laser activates again, i.e. 128 = no cooling. For fused silica 128 or greater is recommended.

#### PULL

Range 0-255 (0 recommended); Generally not required for fused silica.

## Commercial Frit

Commercial frits are available, e.g. NewObjective's IntegraFrit™ see table 1 below for part numbers

OR make your own frit:


## Sintered silica frit

### Reagents and Materials (see Table 1 below)

- Fused silica tubing 360 µm OD by 50, 75, 100, 150, 200 µm ID
- LiChrosorb Si60 5µm underivatized silica (e.g. FisherScientific # M93881)
- Fused silica cutter
- Micro-torch (e.g. small torch for oxygen/propane (Kingsley North Inc.) or Microflame torch)


### Procedure

1. Cut desired length of fused silica
2. Tamp one end of the fused silica into the Lichrosorb Si60, 30 times.
3. Pass this end quickly through the micro torch flame twice to sinter the frit. It should look similar to this frit.

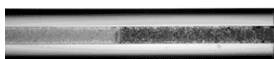


(Image from New Objective's Website)

4. Pack the new frit at 1000 psi and watch to make sure the frit holds

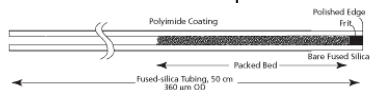

(Image from New Objective's Website)

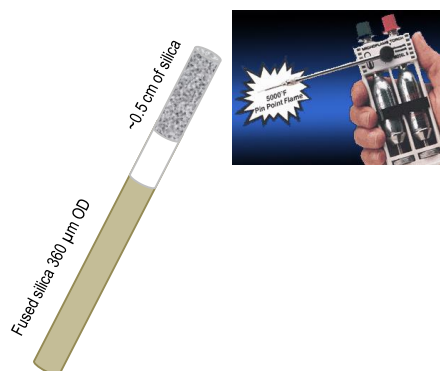

Note the frits can be reused multiple times.

We use the HPLC pump to remove the packing material. Just reverse the fritted fused silica and start the pump at the necessary flow rate, observe the open end of the fused silica until beads flow out and droplet is clear (all the beads are eluted). Check under a microscope to make sure the fritted fused silica is empty.

## Kasil frit


Procedure is based on a method described by Meiring, H.D. et al., J.Sep Sci., 2002, 25, 557-568.

### Reagents and Materials (see Table 1)

- Fused silica tubing 360  $\mu\text{m}$  OD by 50, 75, 100, 150, 200  $\mu\text{m}$  ID
- Fused silica cutter
- KASIL 1 Potassium Silicate Solution (29.1%) or Potassium Silicate, Aqueous Soln, 29.8 Deg. Baume, Pfaltz & Bauer No.:P23830-200g (Fisher part # 50-828-816)
- Formamide, 99.5%, for analysis, ACROS Organics (Fisher part # AC205821000)
- Block Heater or Oven 80-90  $^{\circ}\text{C}$
- Microscope

### Procedure


- Cut fused silica tubing lengths to ~25cm (or longer if needed).
- Use a 1.5 ml Eppendorf tube to make up the frit material this has to be made up and used quickly:
  - Add 200 $\mu\text{L}$  KASIL 1 Potassium Silicate Solution (29.1%) first
  - Add 50 $\mu\text{L}$  Formamide (you can try 100 ul but it may set too fast) and vortex well for a few seconds.
- Or
  - 170 $\mu\text{L}$  Potassium Silicate, Aqueous Soln, 29.8 Deg
  - Add 30 $\mu\text{L}$  Formamide (you can try 100 ul but it may set too fast) and vortex well for a few seconds.
- Spin at 10,000 rpm for 2 min
- Working quickly, place each capillary tube (or as many as you can hold flat between your fingers) into the KASIL/Formamide solution, but not into the pelleted precipitate.
- Capillary action will quickly move the solution up the tube so you only need them in for a brief time. Fill tube with ~1- 2cm of solution.
- Wipe off the outside of the tubing with a kimwipe (all of the bundle can be done at once if it was held flat)
- Inspect under the microscope to see that they are filled to ~1-2 cm to be able to later cut them down to a final length of ~2mm. You want the material to be solid, not broken up as it filled the capillary, though you only care about the portion that will be left.
- Place capillary tubes under a heating block at 90 $^{\circ}\text{C}$ . Leave overnight.  
OR
- Curl the bundle into a ring, but do not bend the tubing at the filled end, as this may fragment the frit material and weaken it. Heat in an oven at 80-85  $^{\circ}\text{C}$ , we normally do about 6 hours this may be overkill.
- Remove tubes from heating block or oven and check each for consistency under a microscope.
- Trim with a tubing cutter before use. About 2 mm gives a good amount of backpressure to the column.
- Using the bomb, pack frit with your packing of choice.

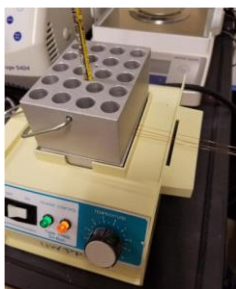

## Column Packing

### Reagents and Materials (see Table 1)

- 80% ethanol (50-100% ethanol, methanol, isopropanol or acetonitrile have successfully been used)
- Packing material: derivatized beads, e.g. Reprosil-Pur 120 C18-AQ, 5  $\mu\text{m}$  (ESI source solutions)
- Fused silica tubing 360  $\mu\text{m}$  OD by 50, 75, 100, 150, 200  $\mu\text{m}$  ID (Polymicro Technologies)
- Column end fittings, frits, and tubing sleeves (optional)
- 1.5ml Eppendorf vials
- Small stir bar
- Microscope and/or light source

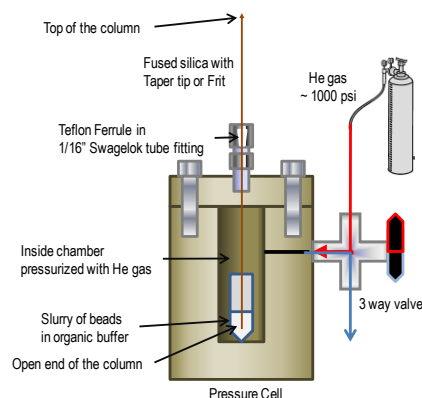

### Procedure

1. Cut a section of silica tubing to desired length and pull a tapered tip or make a frit as described above. The other end of the column will remain "open".
2. To make slurry poor (don't use a spatula to avoid contamination of the stock) a little bit of the packing material (~50  $\mu\text{l}$  dry volume) into an Eppendorf vial and add about 1ml of 80% ethanol. Vortex to mix thoroughly, add a small stir bar in the bottom of the vial. Then cut off the cap with a razor blade and place the vial into the pressure cell chamber. Turn stir plate below cell on. (Note the use of a stirrer is optional)
3. Put the top of the pressure cell in place and tighten the bolts.
4. Thread the "open" end of the capillary into the top of the pressure cell – through the nut, Teflon ferrule and top of the cell. The end of the capillary inside the pressure cell should sit a few millimeters above the bottom of the vial. This is best achieved by pushing the tubing down until you feel it hit the bottom, and then backing it off a bit. When you have it to the correct height, tighten the nut on top of the cell tightly with your fingers. Do not use a wrench for this; it will destroy the Teflon ferrule.
5. Set the high pressure regulator on the He gas tank to about 1000 psi.
6. Slowly increase the pressure in the cell by turning the valve 180° clockwise. You should see a liquid droplet forming at the top of the column.
7. Use a microscope or a light source to observe the beads packing and note the flow of the liquid out the top of the fused silica. The packing process may take 20 minutes-1hour, or longer depending on the inner diameter of the tubing being packed and the length of the final column.
8. Once the desired column length is packed slowly release the pressure from the pressure cell by turning the valve 180° counter clockwise. Remove the column.
9. Connect the column to your LC system and equilibrate with the appropriate buffer. Note the pressure cell can be used to equilibrate the column as well, simply replace the slurry with the appropriate buffer and repeat the steps above.
10. Once the column is equilibrated remove the column and cut it to the exact desired length. Note when using a trap column it is best to avoid any dead volume at the end of the chromatographic column by cutting through the packed part towards the end of the column.

### Trouble shooting

- **Column is not packing, but liquid droplet is forming at the top of the column (most common case)**  
Gently tap the column until beads flow again (may take a minute or so)  
Packing material could be settling: use stir bar to mix slurry or release pressure and vortex slurry again.  
Open end of capillary is clogged, release pressure and check (cut off the end if necessary).
- **No flow is coming out the top during packing**  
Capillary is probably not far enough into the vial, or is clogged. Release pressure and check.
- **slurry is sliding down capillary when pressure is released**  
This is often observed when the pressure of the cell or from the HPLC is released or if the column was not equilibrated long enough.

## Reagents and Materials

Table 1: Reagents and Materials

| Part                    | Description  | Vendor   | part #   |
|-------------------------|--|--|--|
| <b>Torch</b>            | Torch Oxygen/Propane Kit, includes #4 tip, 6' hoses and regulators for disposable propane and oxygen   | <a href="#">Kingsley North Inc.</a>                      | <a href="#">6-1386</a>                           |
|                         | Tip: #4 Replacement Tip (my favorite tip size)   | <a href="#">Kingsley North Inc.</a>                      | <a href="#">6-1391</a>                           |
|                         | Bermzomatic 1.1 Cubic Foot Oxygen Cyli (189-OX9)   | Amazon   | 189-OX9  |
|                         | Propane Burner Replacement Cylinder  | <a href="#">ThermoFisher</a>                             | S41897B  |
|                         | MicroFlame No.4200 - Standard Torch Kit 5000 °F(hard to find)  | <a href="#">Azure Moon Trading Corp</a>                  | <a href="#">No.4200</a>                          |
|                         | Cole Parmer Micro Torch Kit 2500 °F  | <a href="#">ThermoFisher</a>                             | NC9764498  |
| <b>Laser Puller</b>     | P2000 Fiber Puller   | <a href="#">Sutter Instruments</a>                       | <a href="#">P-2000/F</a>                         |
| <b>Packing stations</b> | Homemade pressure cell click on pressure cell on our tools page<br><a href="https://proteomicsresource.washington.edu/protocols05/pressurecell.php">https://proteomicsresource.washington.edu/protocols05/pressurecell.php</a> |  |  |
|                         | Pressure Cell Kit, includes the cell, fittings, manual, tool kit and ferrules (5 x 0.4mm)  | <a href="#">Brechtbuehler Inc.</a>                       | <a href="#">1 1000 100</a>                       |
|                         | Nanobaume™ High Pressure Capillary Packing Assembly, includes the cell, fittings   | <a href="#">Western Analytical Products</a>              | <a href="#">SP-400</a>                           |
|                         | NextAdvance Pressure Injection Cell is rated for 2500 psi. Includes a three-way valve, hex wrench, 10 Ferrules for typical (~360 um OD) capillaries, a frit kit, and operator manual.  | <a href="#">NextAdvance</a>                              | <a href="#">PC77</a><br><a href="#">PC77-MAG</a> |
| <b>Accessories</b>      | Ferrules (replacement) 100% PTFE, 1/16 in. to 0.4mm  | <a href="#">ThermoFisher</a>                             | NC9656755  |
|                         | Regulator 1-Stage High Pressure 3000 PSI Brass, 0-2000 Del Range, CGA-580  | <a href="#">Airgas</a>                                   | <a href="#">Y11N115H580</a>                      |
|                         | Stainless Steel Pipe Fitting, Hex Coupling, 1/4 in. Female NPT   | <a href="#">Swagelok</a>                                 | <a href="#">SS-4-HCG</a>                         |
|                         | SS Swagelok Tube Fitting, Male Connector, 1/8 in. Tube OD x 1/4 in. Male NPT   | <a href="#">Swagelok</a>                                 | <a href="#">SS-200-1-4</a>                       |
|                         | Alltech* Standard Stainless-Steel Tubing, 316, 1/8X085,50FT  | <a href="#">ThermoFisher</a>                             | AT3010   |
|                         | Fisher Scientific* Stereomaster* Microscopes   | <a href="#">ThermoFisher</a>                             | 12-562-12  |
|                         | Fisherbrand* Lab-Jack*   | <a href="#">ThermoFisher</a>                             | 14-673-52  |
|                         | Barnstead/Thermolyne* Cimarec* Digital Stirrers  | <a href="#">ThermoFisher</a>                             | 11-675-910                                       |
|                         | Fisherbrand* Spinbar* Magnetic Micro Stirring Bars   | <a href="#">ThermoFisher</a>                             | 14-513-64  |
|                         | <b>Scribe</b>  | Chromatography Research Supplies INC Column Scribe 10/PK | <a href="#">ThermoFisher</a>                     |
| <b>Frit/Trap</b>        | Silica underivatized (if a frit is needed) EMD CHEMICALS LICRSB SI60 LCH 5UM 10GM  | <a href="#">ThermoFisher</a>                             | M93881   |
|                         | Self-Pack IntegraFrit Columns  | <a href="#">New Objective</a>                            | <a href="#">IF360-100-50-N-5</a>                 |
|                         | KASIL 1 Potassium Silicate Solution (29.1%)  | PQ Corporation PO Box Valley Forge PA 19482              | call for a sample                                |
|                         | Potassium Silicate, Aqueous Soln, 29.8 Deg. Baume, Pfaltz & Bauer No.:P23830-200g  | Fisher   | 50-828-816                                       |
|                         | Formamide - BioUltra, for molecular biology, ≥99.5% (T)  | Sigma-Aldrich  | 47671-250ML-F                                    |
|                         | Formamide, 99.5%, for analysis, ACROS Organics   | Fisher   | AC205821000                                      |
|                         | Block Heater (or Oven 80-90 °C)  | VWR  | 12621-108  |
|                         | BOEKEL Economy Lab Ovens   | BioExpress   | O-2120-07  |
| <b>Fused Silica</b>     | TSP075375 fused silica, 75 µm ID x 360 µOD   | <a href="#">Polymicro Technologies</a>                   | 2000019  |
|                         | TSP100375 fused silica, 100 µm ID x 360 µOD  | <a href="#">Polymicro Technologies</a>                   | 2000023  |
|                         | TSP050375 fused silica, 50 µm ID x 360 µOD   | <a href="#">Polymicro Technologies</a>                   | 2000017  |
| <b>Packing material</b> | Dr Maisch Reprosil-Pur 120 C18-AQ, 5 um (1 gm)   | <a href="#">ESI Source Solutions</a>                     | r15.aq.0001                                      |
|                         | Dr Maisch Reprosil-Pur 120 C18-AQ 3 um (1 gm)  | <a href="#">ESI Source Solutions</a>                     | r13.aq.0001                                      |
|                         | Column Packing Material Jupiter 4u Proteo 90Å, 1 gr  | <a href="#">Phenomenex</a>                               | 04A-4396   |
|                         | Column Packing Material Jupiter 5 u C18 300Å Bulk Packing, 1 gr  | <a href="#">Phenomenex</a>                               | 04A-4053   |

Other vendors of bulk packing material include: Phenomenex, Alltech, Vydac, Restek, ZirChrom... see Table 2 below



Table 2: Bulk Column Packing Media

| Vendor     | Part #       | Description   | Chemistry | Pore Size | Particle size |
|------------|--------------|---|-----------|-----------|---------------|
| Phenomenex | 04A-4396     | Jupiter 4u Proteo 90A, 1 gr   | C12       | 90A       | 4u            |
| Fisher     | AT16082      | Adsorbosphere C18 10 U 80A 10g  | C18       | 80A       | 10u           |
| Michrom    | PM3/91100/00 | Halo 2.7um 90A C18  | C18       | 90A       | 2.7u          |
| Michrom    | PM3/66100/00 | Magic C18 100A 3u PM3/66100/00  | C18       | 100A      | 3u            |
| Michrom    | PM5/66100/00 | Magic C18 100A 5u PM5/66100/00  | C18       | 100A      | 5u            |
| Fisher     | 06-716-307   | Restek Ultra 5microm Bulk Packing Material Ultra C18 Bulk Packing, pore size: 100Å, carbon load: 20%, endcap: fully endcapped, pH range: 2.5 to 7.5 | C18       | 100A      | 5u            |
| Fisher     | 06-717-238   | Pinnacle II C18 Bulk Packing, 5 microm, 110Å, carbon load: 13%, endcap: fully endcapped, pH range: 2.5 to 10, min order 5gr                         | C18       | 110A      | 5u            |
| Michrom    | PM5/66200/00 | Magic C18 200A 5u PM5/66200/00  | C18       | 200A      | 5u            |
| Phenomenex | 04A-4053     | Jupiter 5 u C18 300A Bulk Packing, 1 gr   | C18       | 300A      | 5u            |
| Fisher     | AT16088      | Adsorbosphere C18, 10g 5u   | C18       |           | 5u            |
| Michrom    | PM3/61100/00 | Magic C18AQ 100A 3u PM3/61100/00  | C18AQ     | 100A      | 3u            |
| Michrom    | PM3/61200/00 | Magic C18AQ 200A 3u PM3/61200/00  | C18AQ     | 200A      | 3u            |
| Michrom    | PM5/61100/00 | Magic C18AQ 100A 5u PM5/61100/00  | C18AQ     | 100A      | 5u            |
| Michrom    | PM5/61200/00 | Magic C18AQ 200A 5u PM5/61200/00  | C18AQ     | 200A      | 5u            |
| Michrom    | PM3/67100/00 | Magic C30 100A 3u PM3/67100/00  | C30       | 100A      | 3u            |
| Michrom    | PM5/67100/00 | Magic C30 100A 5u PM5/67100/00  | C30       | 100A      | 5u            |
| Michrom    | PM3/64100/00 | Magic C4 100A 3u PM3/64100/00   | C4        | 100A      | 3u            |
| Michrom    | PM5/64100/00 | Magic C4 100A 5u PM5/64100/00   | C4        | 100A      | 5u            |
| Michrom    | PM3/90100/00 | Halo 2.7um 90A C8   | C8        | 90A       | 2.7u          |
| Michrom    | PM3/60100/00 | Magic C8 100A 3u PM3/60100/00   | C8        | 100A      | 3u            |
| Michrom    | PM3/60200/00 | Magic C8 200A 3u PM3/60200/00   | C8        | 200A      | 3u            |
| Michrom    | PM5/60100/00 | Magic C8 100A 5u PM5/60100/00   | C8        | 100A      | 5u            |
| Michrom    | PM5/60200/00 | Magic C8 200A 5u PM5/60200/00   | C8        | 200A      | 5u            |
| Fisher     | 06-717-157   | Pinnacle II C8 Bulk Packing, 5 microm   | C8        |           | 5u            |
| Michrom    | PM5/62100/00 | Magic Cyano 100A 5u PM5/62100/00  | Cyano     | 100A      | 5u            |
| Michrom    | PM5/63100/00 | Magic Phenyl 100A 5u PM5/63100/00   | Phenyl    | 100A      | 5u            |
| Fisher     | M93881       | EMD CHEMICALS LICRSB S160 LCH 5UM 10GM  | si        | 60A       | 5u            |